

Model Number	3W: SC609 4W: SC709
Weight Capacity	300 lbs.
Maximum Speed ¹	Up to 5.3 mph
Ground Clearance ³	3W: 1.8" at transaxle, 1.4" at anti-tip, 2.7" at side deck 4W: 1.8" at transaxle, 1.4" at anti-tip, 2.1" at motor bracket bolt
Turning Radius ³	3W: 40.25" 4W: 51.5"
Overall Length ³	3W: 43.5" 4W: 45.5"
Overall Width ³	22.25"
Seat-to-Ground Height	19.25"-20.75"
Seat-to-Deck Height	17.25"-18.75"
Front Tires	3W: 3" x 8" solid 4W: 3" x 9" solid
Rear Tires	3" x 9" solid
Suspension	N/A
Range Per Charge ^{1,2}	Up to 13 miles
Total Weight Without Batteries ⁵	3W: 117 lbs. 4W: 131 lbs.
Heaviest Piece When Disassembled ³	3W: 48 lbs. (front section) 4W: 59 lbs. (front section)
Standard Seating ³	Type: Compact, foldable Weight: 28 lbs. Material: Black Vinyl Dimensions: 18" width, 17" depth
Drive System	Rear-wheel drive, sealed transaxle with 24-volt DC motor
Dual Braking System	Electronic, regenerative and electromechanical
Battery Requirements ^{4,6}	Size: (2) U1 standard, Weight: 24 lbs. each
Battery Charger	3.5-amp, off-board
Warranty	Frame: lifetime limited Drivetrain: 2-year limited Electronics: 2-year limited Battery: 6-months

HCPCS Code: K0800 (Group 1, Standard)

- 1 Varies with user weight, terrain type, battery amp-hour (AH), battery charge, battery condition, and tire condition. This specification can be subject to a variance of (+ or -) 10%.
- 2 Tested in accordance with ANSI/RESNA, MC Vol2, Section 4 & ISO 7176-4 standards. Results derived from theoretical calculations based on battery specifications and drive system performance. Test conducted at maximum weight capacity.
- 3 Due to manufacturing tolerances and continual product improvement, this specification can be subject to a variance of (+ or -) 3%.
- 4 Battery weight may vary based on battery manufacturer.
- 5 Includes standard seat weight.
- 6 AGM or gel cell type required.



(US) 800-800-8586 • (Canada) 888-570-1113

pridemobility.com |    

HCPCS codes provided should not be considered as legal advice and do not guarantee reimbursement. DME providers are responsible for determining appropriate billing codes when submitting for insurance reimbursement.

NOTE: All specifications subject to change without notice. The information contained herein is correct at the time of publication; we reserve the right to alter specifications without prior notice. Speed and range vary with user weight, terrain type, battery charge, battery condition and tire pressure. Due to manufacturing tolerances and continual product improvement, length, width, turning radius and ground clearance can be subject to variance of + or - 3%. Special order seats or non-standard dimensions will not be fold-flat and depth-adjustable.

VIC9_111919



VICTORY[®] 9

**COMFORT AND PERFORMANCE
FOR YOUR ACTIVE LIFESTYLE**

LIVE YOUR BEST[®]
pridemobility.com |    


STYLE, COMFORT, PERFORMANCE

18" x 17" Pinchless hinge, black vinyl seat with 8-point swivel and sliders for increased comfort



Standard front basket

Easy-to-grip tiller adjustment knob

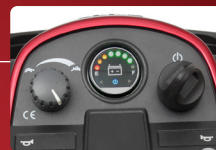


LED battery meter

Long-lasting, maximum-intensity LED headlight offers optimal pathway illumination

Exclusive stylish, lightweight, non-scutting, black, low-profile tires



Tiller console

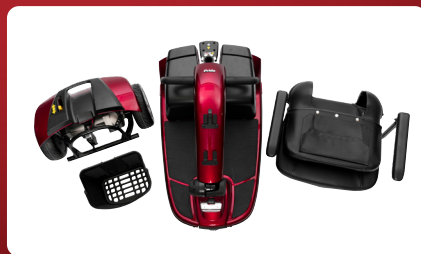

Feather-touch disassembly

The Victory® 9 offers a perfect balance of style, comfort and performance. This scooter has a 300-lb. weight capacity and drives at speeds up to 5.3 mph. With nine-inch solid tires and lots of cool features, the Victory 9 is a great choice!

- Exclusive one-hand feather touch disassembly
- Charger port on the tiller lets you conveniently charge your scooter
- Delta tiller with wraparound handles lets you operate the scooter with one hand and rest your wrist
- (1) standard rearview mirror

Available Colors



Candy Apple Red

Optional Seating

Deluxe contour, high back, black vinyl with height adjustable headrest, flip-up height and width adjustable armrests

- 18" x 18" • 18" x 20"

Accessories

- Cane clip attachment
- Cane/crutch holder (single, double)
- Cell phone holder
- Cup holder
- Forearm crutch holder
- Oxygen tank holder
- Quad cane holder
- Rear basket
- Additional rearview mirror
- Saddlebag
- Safety flag
- Walker holder
- Weather cover
- XLR USB charger
- Lap belts: 50", 60" and 70"

Military Patches

